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1 Safety Information

Warnings

Ensure you read this user guide before installing and using the product.

Babysense is a high quality, high sensitivity electronic device which must be handled with care:

- Do not handle roughly, drop or throw.
- Follow the installation instruction carefully to ensure that it is installed securely and correctly.
- Babysense is not a toy. Do not allow children of any age to play with it.
- Do not put any additional baby or pet in the bed. Babysense cannot be used for two infants in one bed.
- You must undertake the performance test daily of the Babysense monitor (as set out on page 14).
- Do not operate vaporizer or ventilator directly towards the bed.
- To avoid possible danger from sensor cords, ensure that excess cords are pulled back under the mattress, secured with the included cable tie and are not visible or reachable.
- IMPORTANT: It is strongly advised that every caregiver learn how to give first aid.
- Install 4 new and identical size AA 1.5 volt batteries, preferably alkaline.

Please note

"Babysense" monitors and alerts, but cannot prevent cessation of breathing. The responsibility for applying first aid rests on caretakers. The "Babysense" senses breathing movements and motion only, and cannot monitor other occurrences. The purpose of the device is to alert when stopping or slowing of breathing and movement occur.

"Babysense" is a unique safety system which has helped save infants, giving parents (or caretaker) assurance, by enabling intervention in the event of irregular breathing movement or distress, before it becomes critical. There are a number of cautionary steps which are internationally recommended for SAFE SLEEPING:

- Always lay baby to sleep on his or her back.
- Refrain from smoking near or around baby.
- Avoid over-heating of baby's surrounding.
- Do not leave any objects near your baby that can cause suffocation.

Babysense is an addition to these precautions.

"Babysense" is a home breathing movement monitor for basically healthy infants up to one year of age, and is not to be relied on for use with high-risk infants requiring constant medical supervision.

Principle Of Operation

"Babysense" continually detects baby's motion and breathing movements, and will set off a sound and visual alarm if movement stops for 20 seconds or if the breathing rate slows to below 10 breaths per minute. The system's non-touch sensor pads under the mattress monitor baby's movement without restraining him.

Checking the package contents

- Control Unit
- Control Unit Holder
- Optional hanging hook
- 2 Sensor pads
- 2 Cable Ties
- User guide

Please note that the batteries are not supplied.

Disposal

In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to laws on the disposal of **electronic and electrical devices**, owners are obliged to dispose of old electronic and electrical devices in a separate waste container.

The adjacent symbol indicates that the device must not be disposed of in normal domestic waste!

Batteries represent a hazard to both health and the environment! Never open, damage or swallow batteries, nor allow them to get into the environment.

They may contain toxic and ecologically hazardous heavy metals. You are legally obliged to dispose of batteries at the point of sale or at appropriate collection points provided by your local public waste authorities. Disposal is free of charge. The adjacent symbols indicate that batteries must not be disposed of in normal domestic waste and but at authorised collection points.

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
2 Installation

1. Insert, in accordance with the markings, 4 new and identical size AA 1.5 volt batteries, preferably alkaline. Close cover. **Note:** For replacing the batteries, remove the battery cover by inserting a flat screwdriver into the opening recess (see illustration 7). Please dispose of used batteries in accordance with local waste regulations.
2. Place control unit holder on the base of the crib or bed and in between two bars as shown in illustration 1.

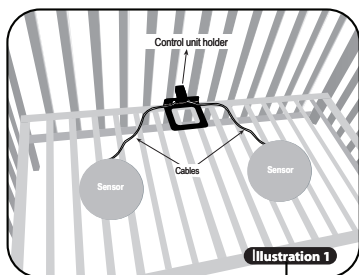
3. Make sure the sensitivity adjuster, located on the back of the sensor pads, is positioned on level 3 (see illustration 2).


4. Sensor pad installation:

- When using for infant not yet able to crawl, place one sensor pad as shown in illustration 3. Make sure the sensor is placed underneath baby's location in crib.
- When baby begins to crawl, and is able to move onto areas not "covered" by a single sensor pad, attach the second sensor and place the two sensors as shown in illustration 3.

PLEASE NOTE: The sensors shall be placed with the marking "THIS SIDE UP" facing upfront.

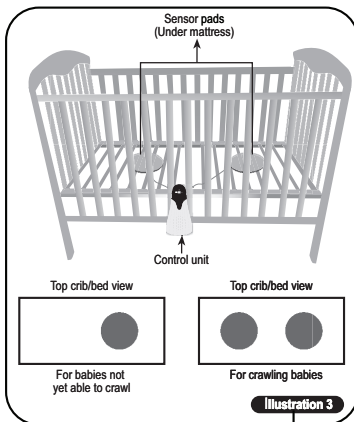
5. Pass the cord(s) of the sensor pad(s) through special openings in the holder. (see illustration 4).

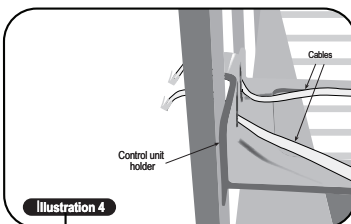
6. Attach end(s) of the cord(s) to the outlets marked "SENSOR" on the back of control unit (see illustration 5 and 6). When using one sensor pad, attach to any "SENSOR" outlet.


PLEASE NOTE: Do not attach any cable to the blocked socket, marked "DO NOT USE".

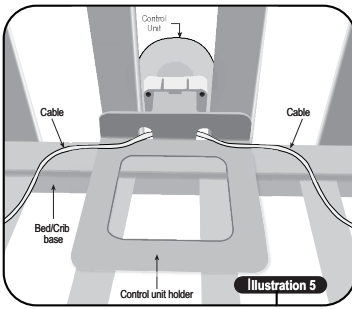7. Should the design of the crib/bed not allow proper installation with the control unit holder, please use the optional hanging hook (see illustration 6) and hang the control unit from the bottom edge of the crib.

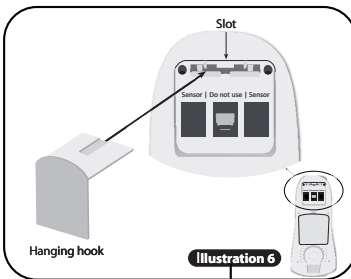8. Press control unit holder or hanging hook firmly into slot in back of control unit (see illustration 5 and 6).

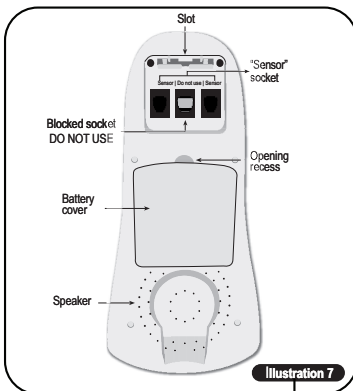
9. Place mattress firmly on top of sensor pads. If the mattress cannot be in contact with the sensor pads, or if the bed base is soft, padded or contains springs, place extra plywood on the bed base, place the sensor pads on the plywood and place the mattress firmly on top.
- IMPORTANT: To avoid possible danger from sensor cords, ensure that excess cords are pulled back under the mattress and are not visible or reachable.**

3 Layout of control unit

See illustration 8 for layout of control unit featuring:

- On/Off button.
- Blue breathing movement indicator light.
- Acoustic and visual red alarm light.
- Low battery red warning light.

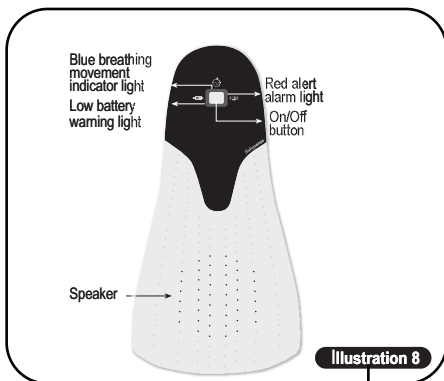4 Operation

"Babysense" is ready for operation and use once the installation is complete. To operate, follow the below instructions:

- Place your baby in crib/bed.
- Turn control unit on by pressing On/Off button.
- Red control lights switch on briefly and a double click sound is heard.
- When breathing and movement are sensed, the blue light flashes.
- Acoustic and visual red alarm light are activated if baby's breathing stops or slows significantly.
- If red low-battery warning light flashes, replace all batteries as soon as possible. If in addition clicking is heard, all batteries must be replaced immediately.
- Turn off by pressing On/Off button before taking your baby out of crib/bed.

Testing Operation

"Babysense" is a highly sensitive device which monitors baby's breathing movements through the mattress during sleep. There are two main causes of malfunction - technical malfunctions and malfunctions due to environmental factors.

Technical malfunction

In order to minimize technical problems, follow installation and operation instructions carefully.

Environmentally caused malfunction

Due to the high sensitivity of "Babysense", it can be influenced by fans, wind, air conditioning, and so forth. Such environmental effects could be picked up by the sensors on a similar level to that of breathing movements, thus preventing or delaying alarm activation.

Conduct this performance test daily

Do not conduct performance test by picking up baby out of bed and waiting for the alarm to sound

- The performance test should be conducted while baby is not in bed.
- Turn device on and take two steps back from bed.
- Ensure that alarm goes off within 20 seconds (due to absence of breathing movements).
- The alarm may sound with or without preliminary ticking.
- Check that alarm can be heard in all locations throughout the home or place of care.
- In case that alarm does not sound within 20 seconds refer to troubleshooting at page .
- Performance test must be conducted daily.

False alarms

Isolated cases of false alarms are possible. Check if caused by faulty connection of the sensor pad(s) to the control unit, connection of only one sensor pad when baby can move off the immediate area, connection of two sensor pads far apart from each other when baby not crawling yet, or failure of the sensor pads to be directly against the mattress.

In case of false alarms, kindly refer to troubleshooting guide. If the problem cannot be solved, please contact your nearest BABYSENSE representative.

In the event of alarm

Act immediately according to the following:

1. Approach the bed and check that baby is breathing. Note pallor (color of skin), muscle tone and movement.
2. If baby does not seem normal, pinch baby's heel and check that there is no cause for choking, such as an object lodged in baby's throat.
When breathing resumes, consult a doctor.
3. If breathing does not resume, apply first aid and seek medical assistance immediately.

Caring for your babysense

"Babysense" requires no special treatment beyond replacing batteries when the battery indicator light goes on. The sensor pads and control unit may be wiped with a damp cloth.

Batteries should be removed and the device stored in a dry shaded place when not in use.

 **It is recommended to replace sensors every two years.**

5 Appendix

Troubleshooting Guide

Should problems arise with the telephone, please refer to the following information first.

In the case of technical problems, contact our Hotline. Please refer to our website

www.audioline.de for the telephone number.

In the case of claims under the terms of guarantee, contact your sales outlet. There is a 2 year period of guarantee.

Symptom - Possible Cause - Solution

Lights on Control Unit don't illuminate when turned ON

- Batteries are not inserted correctly or may be flat.
 - 🔧 Check batteries are inserted in the Control Unit observing the correct polarity.
 - 🔧 Control Unit requires 4 X AA 1.5V Alkaline batteries (not supplied).

Low battery voltage warning (red light) on the Control Unit

- Batteries low.
 - 🔧 Replace all 4 batteries.

False alarm up to 1/2 minute from activating

- Sensor pad is out of order.
- Cable is connected to monitor's blocked socket.
 - 🔧 Connect one sensor pad lay the baby on the bed for 1 minute. If the alarm does not activate, the sensor pad is in order. If the alarm activates, disconnect and place the sensor pad aside.
 - 🔧 Connect the other sensor pad and repeat steps above.
 - 🔧 Do not connect any cables to the blocked socket (only to "SENSOR" socket).

False alarm after 1/2 minute from activating

- The mattress is not touching the sensor pads.
- The mattress is raised (there is an object near the sensor pads).
 - 🔧 Make sure that the mattress and the sensor pads are touching.
 - 🔧 Make sure that there are no objects next to the sensor pads.
- Placing two sensors next to each other (or using one sensor) when the baby is of crawling age or using two sensor pads far apart from each other when baby is young and is not crawling yet.
 - 🔧 Connect according to instructions.
 - 🔧 It is important to use only one sensor pad when the baby is young and is not crawling yet.
 - 🔧 It is important to make sure that two sensor pads are used when the baby is moving around.
 - 🔧 Make sure the sensors are placed underneath baby's location in bed/crib.
- The mattress does not enable transfer of baby movements. (a rigid mattress, a mattress with springs or a mattress with air pockets).
 - 🔧 Place extra plywood under the sensor pad so that the sensor pad and the mattress will be tightly pressed together.

Delayed alarms while taking the baby out of the bed

- Performance test is not correct.
 - 🔧 Do not conduct performance test by picking up the baby out of bed and waiting for the alarm to sound. Kindly conduct performance test according to guidance in page 14.

The alarm does not go off within 20 seconds when performing the correct performance test

- Environmental influences detected by the sensor pads, such as wind, fan, air conditioner, mobile, vibrations.

Perform a test which will detect the source of the movement that the device detects. For example, sources of wind such as: air conditioner, an open window or a fan; movements of a baby mobile (which can vibrate even if it's not turned on) walking around the bed area or an intercom placed near the bed with intense transmission levels. Solutions: Eliminate sources of wind, remove the mobile etc. After neutralizing all environmental influences, perform the test again. If even after neutralizing all of the environmental influences, the alarm does not go off within 20 seconds, perform the following:

1. In the back of the sensor pads (see illustration 2) reduce the sensitivity level to position 2 (on both sensor pads) and perform again the Performance Test. If the alarm sounds within 20 seconds, than keep this sensitivity level.
2. If the alarm does not go off in time, than reduce the sensitivity level to position 1 and perform again the Performance Test. If the alarm sounds within 20 seconds, than keep this sensitivity level.
3. If the alarm does not go off in time at this level, than it is suggested using only one sensor pad with the sensitivity level adjusted on position 1.
4. **PAY ATTENTION:** Should false alarms sound as a result of reduced sensitivity, than increase the level of sensitivity and perform the Performance Test.

Technical Specifications

Powered by 4 size AA, 1.5 volt batteries.

Power consumption approx. 1 mA.

Battery life approx. 6 months using alkaline batteries.

Alarm volume typically 85dB SPL at distance of one meter.

Size of sensor pad is ø 216 mm.

Declaration of Conformity

This device fulfils the requirements stipulated in the EU directives.

BABYSENSE Infant Movement Respiratory Monitor (Model CU-100/2 revD.) is in conformity with the essential requirements and provisions of Annex VI and Annex III of the Medical Devices Directive 93/42/EEC as amended by 98/79/EC, 2000/70/EC, 2001/104/EC and 2007/47/EC.

The product is classified as Class IIb medical device, in accordance with Rule 10, Annex IX of the Directive.

Notified Body:

LNE No. 0459

1 rue Gaston Boissier

75724 PARIS Cedex 15, France

Authorised Representative:

Obelis S.A.

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1030 Brussels, BELGIUM

Conformity with the above mentioned directives is confirmed by the CE symbol on the device.

To view the complete Declaration of Conformity, please refer to the free download available on our web site **www.audioline.de**.

Guarantee

AUDIOLINE equipment is produced and tested according to the latest production methods. The implementation of carefully chosen materials and highly developed technologies ensure trouble-free functioning and a long service life. The terms of guarantee do not apply where the cause of equipment malfunction is the fault of the telephone network operator or any interposed private branch extension system. The terms of guarantee do not apply to the batteries or power packs used in the products. The period of guarantee is 24 months from the date of purchase.

All deficiencies related to material or manufacturing errors within the period of guarantee will be redressed free of charge. Rights to claims under the terms of guarantee are annulled following tampering by the purchaser or third parties. Damage caused as the result of improper handling or operation, normal wear and tear, incorrect positioning or storage, improper connection or installation or Acts of God and other external influences are excluded from the terms of guarantee. In the case of complaints, we reserve the right to repair defect parts, replace them or replace the entire device. Replaced parts or devices become our property. Rights to compensation in the case of damage are excluded where there is no evidence of intent or gross negligence by the manufacturer.

If your device shows signs of defect during the period of guarantee, please return it to the sales outlet in which you purchased the AUDIOLINE device together with the purchase receipt. All claims under the terms of guarantee in accordance with this agreement can only be asserted at the sales outlet.

No claims under the terms of guarantee can be asserted after a period of two years from the date of purchase and hand-over of the product.

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

Operating Elements	2
Starting Up	3
Checking the package contents	3
Connecting and aligning the transmitter	3
Operating with batteries (not supplied)	3
Operating with a power adapter plug	3
Connecting the receiver	4
Operating with a power adapter plug	4
Operation with rechargeable batteries	4
Mounting/dismounting the belt clip	4
Optimization of range	4
Operating the Transmitter	5
Switching the transmitter on/off	5
Connecting and aligning the transmitter	5
Night vision	5
Operating the Receiver	6
Switching the receiver on/off	6
Setting the volume	6
Setting the brightness	6
Temperature option (setting the display unit)	6
Temperature option (setting the range)	6
Switching the night light on/off	7
Switching the vibration alarm on/off	7
Adjusting the VOX sensitivity	7
Video/Audio switching	7
Setting the zoom	8
Selecting a lullaby	8
Intercom function	8
Range	8
Registration	8
Multi camera mode	9
External shaker option	9
Appendix	10
Safety Information	10
Intended use	10
Installation location and safety information	10
Power adapter plug	10
Rechargeable batteries	10
Power failure	10
Medical equipment	11
Disposal	11
Helpline	11
Troubleshooting guide - Problems and solutions	11
Declaration of conformity	12
Maintenance	12
Technical data	13
Guarantee	13

Operating Elements

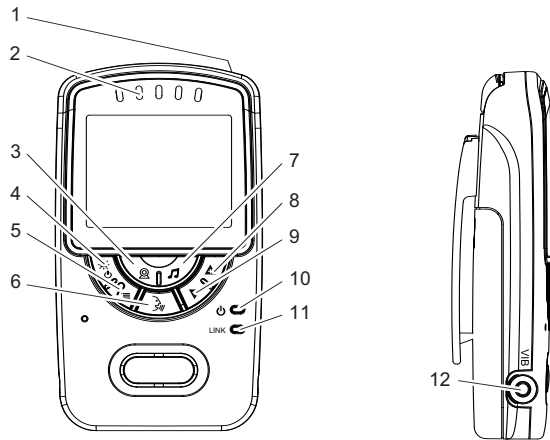
Monitor unit (receiver)

- | | | | |
|---|---------------------------------------|----|---|
| 1 | Antenna (fold-out) | 7 | Lullaby (on/off / selection) |
| 2 | Sound level LEDs | 8 | Increase volume / menu options selection |
| 3 | Camera selection / Registration | 9 | Reduce volume / menu options selection |
| 4 | On/Off switch / Video/audio switching | 10 | Power LED |
| 5 | Menu button / Zoom on/off | 11 | Link LED |
| 6 | Talk button | 12 | Socket for external vibration alert device
(e. g. Amplicomms PTV 100, refer to page 9) |

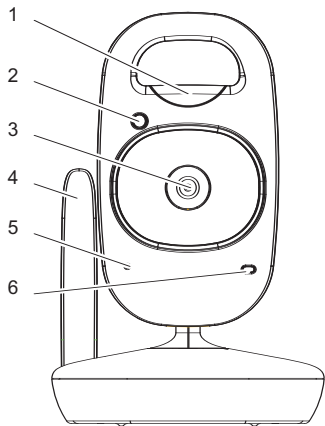Camera unit (transmitter)

- | | | | |
|---|--------------|---|-------------------------------|
| 1 | Night light | 6 | Power LED |
| 2 | Light sensor | 7 | ON/OFF switch |
| 3 | Camera lens | 8 | Socket for power adapter plug |
| 4 | Antenna | 9 | Temperature sensor |
| 5 | Microphone | | |

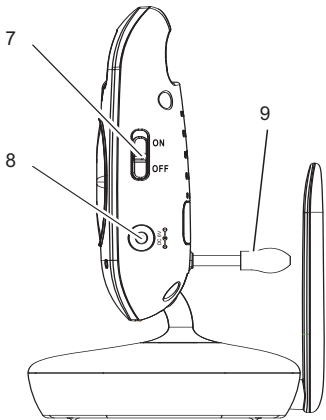Starting Up

⚠ Important: please make sure you've read the safety information on page 10 before you set up your system. The data transmitted (image and/or sound) are digitally coded.

Checking the package contents


The package contains:	V130	V132
• One camera unit (transmitter)	1	2
• One monitor unit (receiver) with belt clip	1	1
• One rechargeable battery pack for the receiver	1	1
• One power adapter plug for the transmitter	1	2
• One charger with power adapter plug for the receiver	1	1
• One operating manual	1	1

Connecting and aligning the transmitter

⚠ Attention: Position the transmitter with a minimum distance of 1 m to other electronic devices, otherwise there is a risk of mutual disturbance. Point the upper section of the transmitter towards the baby or object you want to monitor. Check that the baby or object is suitably displayed on the receiver.

Operating with batteries (not supplied)

Use a small screwdriver to open the transmitter's battery compartment on the underside. Insert the batteries, paying attention to correct polarity. Refer to the diagrams in the battery compartment. Close the battery compartment again.

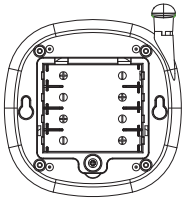
The transmitter is operated by means of 4 AAA Alkaline batteries (not contained in the material supplied). The transmitter cannot be charged. Insert the batteries before connecting the power adapter plug.

In the case of a low battery charge, the Power LED flashes slowly. The batteries must be replaced.

Operating with a power adapter plug

Connect the transmitter as illustrated in the diagram. For safety reasons, only use the power adapter plug supplied.

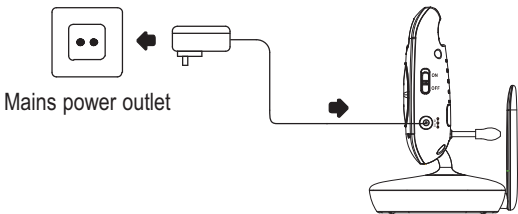
Starting Up

Connecting the receiver

⚠ Attention: Position the receiver with a minimum distance of 1 m to other electronic devices, otherwise there is a risk of mutual disturbance.

Operating with a power adapter plug

Connect the receiver as illustrated in the diagram. For safety reasons, only use the power adapter plug supplied.

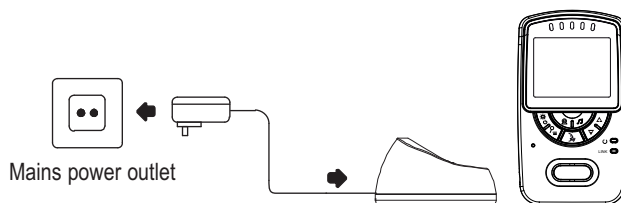Operation with rechargeable batteries

1. Open the battery compartment on the rear side of the receiver.
2. Insert the battery pack, and connect it with the socket. The connector can only be fitted one way.
3. Close the battery compartment again.


The receiver is powered by a Lithium-ion battery pack of 3.7 V, 1500 mAh.

The rechargeable battery pack can be recharged directly in the receiver when it is inserted in the charging station. The Power LED lights up when the receiver is being charged.

In the case of a low battery charge, the Power LED flashes slowly. The battery pack must then be recharged.


Mounting/dismounting the belt clip

Prize the latch (see figure) with a suitable tool (such as a small screwdriver) a little bit away from the chassis and move the belt clip down simultaneously. To attach the belt clip slide it into the guide towards the top until it snaps in.

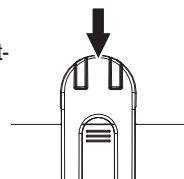Optimization of range

Baby unit and parent unit come with fold-out antennas. If necessary, optimum range can be achieved if both the Baby unit and Parent unit antennas are extended vertically.

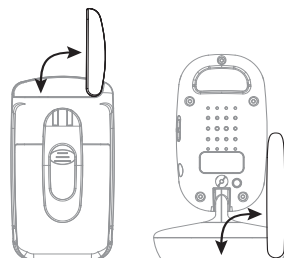Operating the Transmitter

Check the transmitter and receiver are working properly each time before the equipment is put into operation.

Switching the transmitter on/off

Position the transmitter in a good position for the monitoring tasks required and align the camera lens to the area to be monitored.

Slide the switch to ON. The Power LED lights up.
The transmitter is now ready to operate.

To switch the unit off, slide the switch to OFF.

In the case of a low battery voltage, the Power LED flashes slowly.
In this event replace the 4 x AAA batteries with fresh alkaline cells.

Connecting and aligning the transmitter

 **Attention:** Position the transmitter with a minimum distance of 1 m to other electronic devices, otherwise there is a risk of mutual disturbance. Point the upper section of the transmitter towards the baby or object you want to monitor. Check that the baby or object is suitably displayed on the receiver.

Night vision

The transmitter is equipped with 6 infrared diodes. This enables the transmission of pictures from areas which are not illuminated. The night vision function is activated automatically when the light conditions make it necessary.

When the night-time vision function is activated, the screen display on the receiver turns to black and white.

Operating the Receiver

Operating the Receiver

Check the transmitter and receiver are working properly each time before the equipment is put into operation.

Switching the receiver on/off

Press and hold this button until the logo appears first then the image being taken by the transmitter camera.

Press and hold the button for approx. 1 second to switch the unit off.

In the case of a low battery charge, the Power LED flashes slowly. Place the receiver in the charger to recharge the battery pack.

Setting the volume

Use these buttons to adjust the volume to one of five settings or "Off".

Setting the brightness

Press the button once.

This icon appears in the display.

Use these buttons to adjust the brightness to one of five settings.

Temperature option (setting the display unit)

In the lower right hand corner the temperature in the baby's room is shown.

IMPORTANT INFORMATION The ideal temperature for a baby to sleep is 16° - 18 °C. Please set or adjust the room temperature accordingly.

Press the button twice.

This icon appears in the display.

Use these buttons to set temperature measurement to °C or °F.

Temperature option (setting the range)

You can adjust the upper and lower limit for the temperature monitor. An alert will occur if the temperature becomes outside the preset temperature range.

Press the button three times (lower limit) or four times (upper limit).

This icon appears with the current setting in the display (Lo = lower limit / Hi = upper limit).

Use these buttons to adjust the limits (lower limit Off, 12, 15, 18 / upper limit: 20, 23, 26, Off). When set to "Off" the temperature monitoring is deactivated.

Switching the night light on/off

Press the button five times.

This icon appears in the display.

Use these buttons to switch the night light on the transmitter on or off.

Switching the vibration alarm on/off

Press the button six times.

This icon appears in the display.

Use these buttons to switch the vibration alarm on or off.
When the vibration alarm is activated, the receiver vibrates when sounds in the room being monitored are so loud that the red Sound Level LEDs are activated.

Adjusting the VOX sensitivity

(VOX: Voice Operated Transmission = Sound detection)

The VOX function (sound detection) causes the system to start transmitting automatically from the transmitter to the receiver immediately a voice or sound is detected.

Press the button seven times.

This icon appears in the display.

Use these buttons to set VOX to OFF or adjust the sensitivity to one of 5 levels (5 = highest sensitivity). When one of the five levels has been set, the VOX icon lights up in the display.

If the receiver is placed on the charger the display will switch off after approx. one minute of inactivity, but will switch back on if triggered by the vox detector in the transmitter or by pressing any button on the receiver.

However if VOX is set to OFF, the display will remain on all the time.

Video/Audio switching

You can manually switch off the display to save battery power. The sound from the baby's room are to be heard permanently and the display will automatically switch on again if triggered by the vox detector.

Press this button briefly.

For switching the display on again press this button again briefly.


Operating the Receiver

Setting the zoom

You can zoom in on the image.

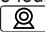Press and hold this button for approx. 2 seconds.

The zoom will be activated. You can move the selection, as long as four arrows are shown in the display (approx. 5 seconds), by pressing the button  to the left,  to the right,  upwards or  downwards.

To deactivate the zoom function press and hold this button again for approx. 2 seconds.

Selecting a lullaby

Switch on: press the button briefly.

This icon appears in the display.

Use these buttons to select one of the 5 lullabies (1-5).

Switch off: press the button briefly.

Intercom function

Press and hold the button.

This icon appears in the display.

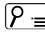
Speak into the microphone to talk with the person near the other unit. When the talking button is pressed, any playing of lullabies will be aborted.

Range


When the distance between the transmitter and receiver is too far, the message "Out of Range" appears in the display and you hear a beep. Move the receiver nearer to the transmitter. If necessary, optimum range can be achieved if both the Baby unit and Parent unit antennas are extended vertically.

Registration

In delivery state the transmitter is pre-registered at the receiver (Watch&Care V132: both transmitters). If necessary or if a another camera is added to the system, registration can also be completed manually:

- 1 Switch the transmitter off (switch set to **OFF**). Remove the batteries of the receiver for a few seconds if necessary.
- 2 Switch the receiver on.
- 3 Press the menu key repeatedly  until [C1] appears in the display.
- 4 Use the   buttons to select the desired camera [C1] ... [C4].
- 5 Press the  button briefly. The Link LED flickers for approx. 10 seconds. During this time, you must perform steps 6 and 7.
- 6 Switch the transmitter on (switch set to ON).

- 7 Press the registration key on the back side of the transmitter. After a short time the video from the transmitter will appear on the screen if the registration was successful.
If the registration failed, "Out of Range" appears on the screen.

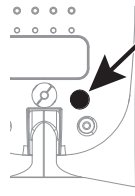Multi camera mode

If you registered several cameras (up to 4) to your system, you can select a single camera or switch to scan mode. In scan mode, the signals from all activated cameras are continuously displayed alternately each for about 10 seconds.

Press this button briefly. The number of the selected camera is shown.

Select another camera or the scan mode by pressing this button repeatedly.

Press this button to confirm your selection.

Selecting the cameras for the scan mode

When performing the selection of the camera (see above), you can set by pressing these buttons for each of the registered cameras, whether it should be assigned to the scan mode (e. g. "C2 x" = not assigned to the scan mode / "C2 ✓" = assigned to the scan mode).

For each of the four cameras you can select the features settings individually.

External shaker option

1. The external shaker (Amplicomms PTV100) is an optional device which is supplied separately to this system. The external shaker can be plugged into the designated output at the side of the display on the parent unit, to enable use as a powerful under-pillow vibrating alert. The external shaker operates when the parent unit built-in vibrator does. Therefore the vibrator option must be enabled within the menu.
2. The parent unit display and sound, is activated when the system detects the baby crying, the sensitivity level of which can be adjusted in the VOX setting. A single activation of the external vibrator/shaker also occurs. The parent unit stays active for 60 seconds following the trigger, during which time any further triggering of the vibrator/shaker will only occur when either of the two sound level LEDs are activated by the incoming sound. Therefore the incoming sound level has to be reasonably loud (at least 75dB@0.5mtrs) in order to trigger the 2nd highest LED. This is set to avoid false triggering the system by one-off low level sounds such as baby sighing for example.

Appendix

Safety Information

Please read this operating instruction manual thoroughly.

Intended use

The surveillance system has been developed to monitor rooms. The transmitter sends acoustic and video signals from its point of installation to the receiver. Any other use is considered unintended use. The use of this equipment does not substitute personal supervision, e.g. of a child. Users are not relieved of their personal liability. Unauthorised modification or reconstruction is not permitted. Under no circumstances open the device or complete any repair work yourself.

Installation location and safety information

- Never lay the transmitter in a baby's/child's bed or playpen.
- Ensure that the transmitter and power cable are out of reach of babies/children, i.e. they are at least one meter away.
- Pay attention that the components are sufficiently well-ventilated when in operation. Do not lay a cushion, towel or such over the devices.
- The devices are not intended for medical purposes. Warning and signal tones from respiratory or heartbeat monitoring equipment are not transmitted.
- The installation location plays a decisive role in ensuring proper operation. Therefore, maintain a distance of at least one meter from other electronic equipment, such as microwave ovens or hi-fi devices, otherwise they could cause mutual interference/faults.
- Prevent excessive exposure to smoke, dust, vibration, chemicals, moisture, heat and direct sunlight. Do not use the equipment in potentially explosive areas.

Power adapter plug

 **Attention:** Only use the power adapter plugs contained in the material supplied since other power adapter plugs could damage the surveillance system. Ensure access to the power adapter plug is not obstructed by furniture or such.

When using the equipment with a power adapter plug, pay attention that the plug and cable are in perfect condition. Kinked or worn cable represents the risk of a fatal accident!

The power adapter plug supplied fulfils the ecodesign requirements of the European Union (Directive 2005/32/EC). This means that, both in an operating state and in an idling state, the power consumption is considerably lower compared to power adapter plugs with an older design.

- Énergiesparend
- Economisant l'énergie
- A risparmio di energia
- Energy efficient

Rechargeable batteries

 **Attention:** Never throw batteries into a fire. Only use batteries of the same type! Pay attention to correct polarity.

Power failure

In the event of a power failure, the devices only operate through battery power. Disconnect the power adapter plugs from the devices in this case.

Medical equipment

 **Attention:** Never use the devices in the vicinity of medical equipment. Effects on such equipment cannot be fully ruled out.

Disposal

In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to laws on the disposal of electronic and electrical devices, owners are obliged to dispose of old electronic and electrical devices in a separate waste container. The adjacent symbol indicates that the device must not be disposed of in normal domestic waste!

Batteries represent a hazard to health and the environment! Never open, damage or swallow batteries or allow them to pollute the environment. They may contain toxic, ecologically hazardous heavy metals. You are legally obliged to dispose of power packs and batteries at the point of sale or in the corresponding containers provided at collection points provided by local public waste authorities. Disposal is free of charge. The adjacent symbols indicate that the batteries must not be disposed of in normal domestic waste and must be brought to collection points provided by local public waste authorities.

Packaging materials must be disposed of according to local regulations.

Helpline

Should problems arise with the device, please refer to the following information first.

In the case of technical problems, contact our Hotline. Please refer to our website www.audioline.de for the telephone number.

In the case of claims under the terms of guarantee, contact your sales outlet.

There is a 2 year period of guarantee.

Troubleshooting guide - Problems and solutions

Equipment does not function

- Transmitter and/or receiver are not switched on.
- The batteries are empty or defective.
- The power adapter plugs are not connected.

No reception

- Transmitter and receiver are not registered.
- Transmitter is too far from object being monitored, e.g. baby.
Observe the minimum distance of one meter.
- The receiver and transmitter are too far apart.
- The VOX sensitivity is set too low.

Reception interference

- The receiver and transmitter are too far apart.
- The battery power is too low.
- Other electronic devices are causing interference in the field of transmission.

Howling or squealing noise (acoustic feedback)

- Increase the distance between receiver and transmitter.
- Reduce the volume on the receiver.

Noise

- The battery power is too low.

Picture but no sound transmission

- Volume setting is too low.

Picture is in black and white

- The night vision function is active.

The battery pack in the receiver is not charged

- The power adapter plug is not connected.

Declaration of conformity

This device fulfils the requirements stipulated in the EU directive:

1999/5 EU directive on radio equipment and telecommunications terminal equipment and the mutual recognition of their conformity.

Conformity with the above mentioned directive is confirmed by the CE symbol on the device.

To view the complete Declaration of Conformity, please refer to the free download available on our website **www.audioline.de**

Maintenance

- Clean the housing surfaces with a soft, fluff-free cloth.
- Never use cleaning agents or solvents.

Technical data

Camera (transmitter)	Socket: 6 V DC Sensor: 1 night light sensor IR diodes: 6
Monitor (receiver)	Screen: 2,4" LC display Socket: 6 V DC
Power supply (receiver)	Power adapter plug, input: 100-240 VAC, 50-60Hz, 150 mA Power adapter plug, output: 6 VDC, 800 mA or Lithium-ion battery pack 3.7 Volt / 1500 mAh
Initial charging time for battery	Up to 16 h
Power supply (transmitter)	Power adapter plug, input: 100-240 VAC, 50-60Hz, 150 mA Power adapter plug, output: 6 VDC, 800 mA or 4 AAA Alkaline batteries
Range	Up to 150 meters under optimum conditions. Range is reduced in buildings through walls, doors etc.
Frequency range	2.4 GHz (digital)
Permissible ambient temperature	0°C to 40°C
Transmission power (TX PA power)	15 ~ 17 dBm

Guarantee

AUDIOLINE equipment is produced and tested according to the latest production methods. Carefully selected materials and highly developed technology ensure perfect functioning and a long service life. The terms of guarantee do not apply to the batteries or power packs used in the products. The period of guarantee is 24 months from the date of purchase.

All deficiencies related to material or manufacturing errors within the period of guarantee will be redressed free of charge. Rights to claims under the terms of guarantee are annulled following tampering by the purchaser or third parties. Damage caused as the result of improper handling or operation, incorrect positioning or storage, improper connection or installation or Acts of God and other external influences are excluded from the terms of guarantee. In the case of complaints, we reserve the right to repair defect parts, replace them or replace the entire device. Replaced parts or devices become our property. Rights to compensation in the case of damage are excluded where there is no evidence of intent or gross negligence by the manufacturer.

If your device does show signs of a defect within the period of guarantee, please contact the sales outlet where you purchased the AUDIOLINE device, producing the purchase receipt as evidence. All claims under the terms of guarantee in accordance with this agreement can only be asserted at the sales outlet. No claims under the terms of guarantee can be asserted after a period of two years from the date of purchase and hand-over of the product.

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AUDIOLINE GmbH